CD40

 Chemdye®
Rev. 15 / 07.2024

Multivariable Process Indicators

Result Reference Guide


Unprocessed


Failures

One or more sterilization parameters were not achieved. The materials cannot be considered sterile and must be processed again. The main failure causes are:

- Inadequate cycle duration
- Inadequate cycle selection
- Overloading / incorrect loading of the sterilizer
- Improper sterilization agent concentration
- Sterilizer malfunction
- Blocking indicator with large metal mass

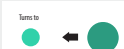Processed and Satisfactory

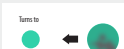Processed and Satisfactory*

***NOTE:** This result applies to uses outside the USA where exposures could be prolonged in time, concentration and temperature beyond those specified in the intended use systems.

NOTE: Reference colors exhibited in printed prospects as well as those displayed in our website and soft copies of other documents, do not necessarily represent the real color as shown in the actual indicators.

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Explanation of Symbols

Product designed for use with Plasma or Vaporized Hydrogen Peroxide sterilization cycles.
CE mark.
Authorized representative in the European Community.

Medical Device.
Batch number.
Expiration Date.
Manufacturer.

Intended Use systems	
Model	Vaporized Hydrogen Peroxide
	Cycles
CD40	Sterrad® 100S Sterilization System Sterrad® NX Sterilization System (Standard and Advanced cycles)
	Sterrad® NX with ALLClear® Technology Sterilization System (Standard and Advanced cycles)
	Sterrad® 100NX Sterilization System (Standard, Flex, Express and Duo cycles)
	Sterrad® 100NX with ALLClear® Technology Sterilization System (Standard, Flex, Express and Duo cycles)
	V-Pro® maX Low Temperature Sterilization System (Non Lumen, Lumen and Flexible cycles)
	V-Pro® maX 2 Low Temperature Sterilization System (Fast Non Lumen, Non Lumen, Lumen and Flexible cycles)

Terragene® and Chemdye® are registered trademarks of Terragene S.A.
STERRAD® is a registered trademark of Advanced Sterilization Products.
V-PRO® is a registered trademark of Steris Corporation.

ATTENTION: If your chemical indicator has this logo printed on it, it means that your product is compatible with exclusive **Trazanto** Automatic Quality Control and Traceability System for cleaning, disinfection and sterilization processes.
This System interprets the indicator's results in an accurate and automatic manner, and stores all the process data through Bionova® Software. Use Trazanto scanner to interpret the results (see instructions for use in the manual of the device).
Do not write on the indicator before scanning.

EN Multi-Variable Chemical Indicator Strip

For Plasma or Vaporized Hydrogen Peroxide sterilization

Indications for Use
United States
Terragene® Chemdye® CD40 is a multivariable chemical indicator intended for monitoring the efficacy of vaporized hydrogen peroxide sterilization processes. The chemical indicator changes from purple to green to indicate that the conditions of the cycle have been met. See Intended Use systems for reference.

Outside the United States
Terragene® Chemdye® CD40 has been designed for monitoring vaporized hydrogen peroxide sterilization processes.

Product Description
Terragene® Chemdye® CD40 (Type 4 according to ISO 11140-1:2014 standard) are single-use chemical indicators that consist of plastic strips printed with indicator ink. These indicators have been designed to monitor plasma or vaporized hydrogen peroxide sterilization processes within loads, ensuring an adequate exposure to the sterilizing agent during the sterilization process.
Innovative highly sensitive purple ink was developed to turn to green or lighter when the indicator is exposed to plasma or vaporized hydrogen peroxide considering the critical process variables which provides additional quality assurance.

Indicators 100 % free of toxic heavy metals.

Precautions
Do not expose Chemdye® CD40 to ethylene oxide (EO), dry heat or any sterilization process other than plasma or vaporized hydrogen peroxide. Do not reuse.

Instructions for Use
1. Place a Chemdye® CD40 within each package, bag or tray that will be sterilized by Plasma or Vaporized Hydrogen Peroxide. In rigid containers, place the indicator in every corner, or at least two diagonally opposite corners, to ensure that the sterilization agent has reached the most inaccessible areas.
2. Perform the sterilization cycle.
3. After the sterilization process has finished, remove the indicators from the sterilizer and analyze the results. For all the chemical indicators exposed to plasma or vaporized hydrogen Peroxide, the ink must turn to the green or lighter reference color when the sterilization parameters are reached. For reference color, please refer to the Result Reference Guide.
NOTE: If any serious incident occurs in relation to the device, it should be reported to Terragene S.A. and the competent authority of the State in which the user is established.

Monitoring frequency
Follow the policies and procedures with the monitoring frequency specified by the professional associations and/or standards corresponding to your country. As recommended practice, and to provide optimal patient safety, Terragene S.A. recommends that each sterilization load be monitored with the appropriate biological indicator.

Storage
Indicators should be stored in a dry place, protected from light at temperatures between 10-30 °C (50-86 °F), 30-80 % relative humidity. Do not get wet. Do not store near sterilizing agents.

Shelf Life
Chemical indicators have an expiration date of 5 years from the date of manufacture when stored at recommended storage conditions. Both the expiration date and batch number of the product are specified on the package label and in each unit of the product. Do not use indicators after their expiration date.
Endpoint Stability Reaction: chemical indicator endpoint shall remain unchanged for a period of not less than 6 months when stored at previously indicated conditions.

Disposal
Discard chemical indicators after use along with plastic waste, according to your country's healthcare and safety regulations.

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