



ADVANCED PIPET SYSTEMS

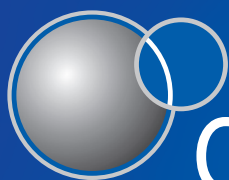

COPAN

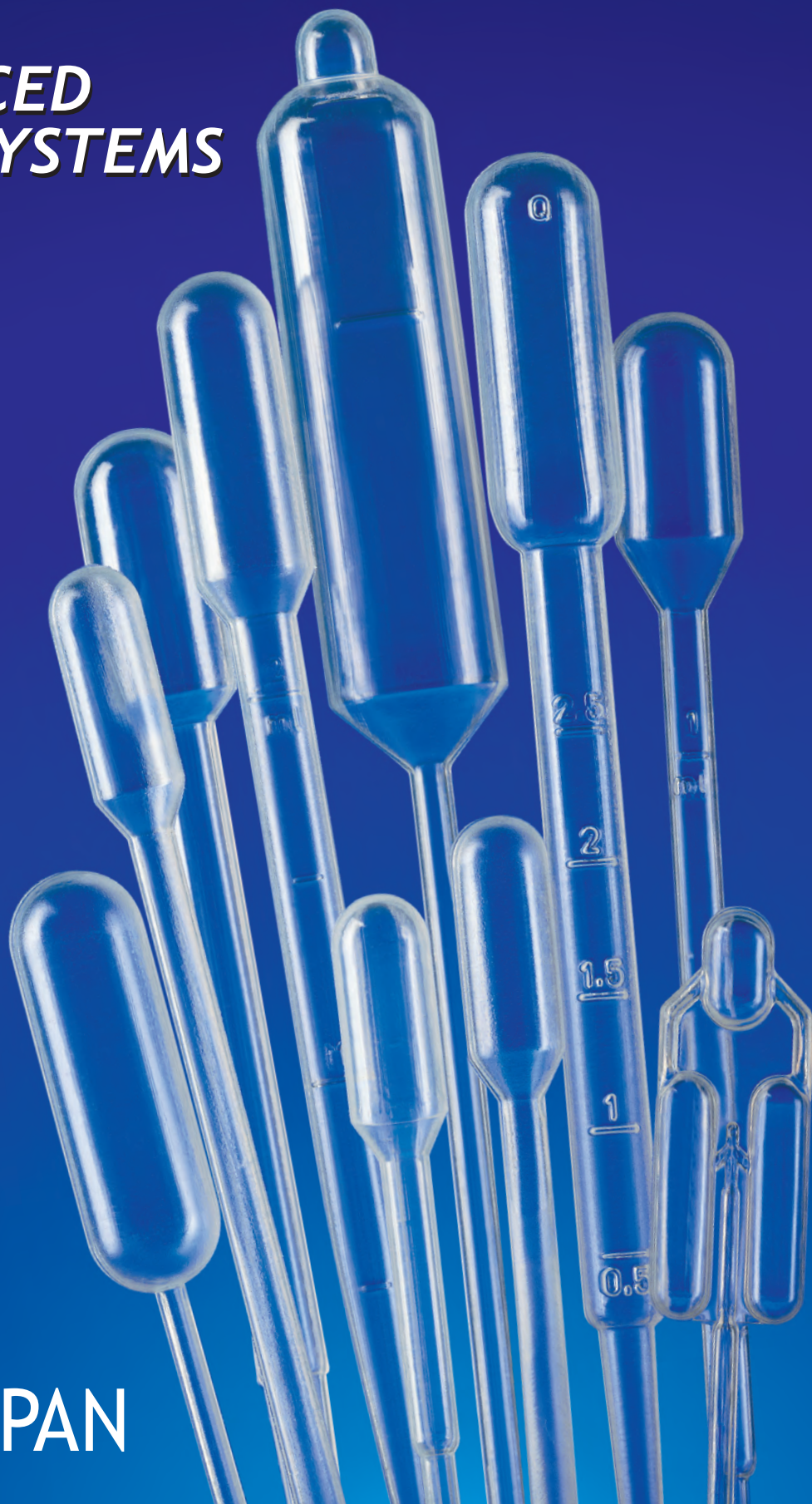

COPAN DISPOSABLE PASTEUR PIPETS

- Ideal for transferring and dispensing liquids safely in all types of laboratories (Blood Bank, Haematology, Bacteriology, Urinalysis, etc.) eliminating the risk of cross contamination.
- Produced in low density polyethylene - they are non-toxic and can be incinerated.
- Can be heat sealed and frozen at very low temperatures, down to -70°C.
- When sterile, packages are in medical grade paper and atoxic film.
- Available in a wide range of standard and special types.


GRADUATED PIPETS



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw	Graduated up to	Grad. intvl.
150 mm	7.8 mm	21	7 ml	3.5 ml	3 ml	0.5 ml
Code	Description	Quantity per box		Quantity per case		
200C	non sterile	500 pieces		5000 (10 boxes of 500 pcs.)		
200CS01	sterile individually wrapped	500 pieces				
200CS20	sterile 20 pieces per pack	25 packs of 20 (500 pcs total)				



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw	Graduated up to	Grad. intvl.
150 mm	5.0 mm	23	5 ml	3.5 ml	1 ml	0.25 ml
Code	Description	Quantity per box		Quantity per case		
201C	non sterile	500 pieces		5000 (10 boxes of 500 pcs.)		
201CS01	sterile individually wrapped	500 pieces				
201CS20	sterile 20 pieces per pack	25 packs of 20 (500 pcs total)				

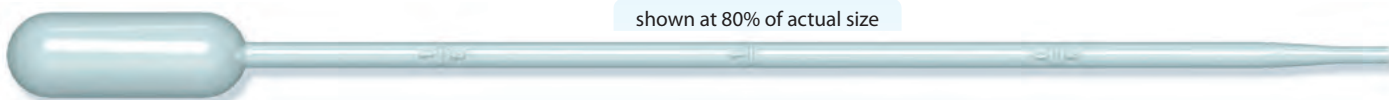

shown at 80% of actual size

Total length	Stem diameter	Drops/ml	Capacity	Bulb draw	Graduated up to	Grad. intvl.
230 mm	5.0 mm	22	6 ml	3.5 ml	1.5 ml	0.25 ml
Code	Description	Quantity per box		Quantity per case		
204C	non sterile	400 pieces		2400 (6 boxes of 400 pcs.)		
204CS01	sterile individually wrapped	500 pieces				
204CS20	sterile 20 pieces per pack	25 packs of 20 (500 pcs total)				



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw	Graduated up to	Grad. intvl.
138 mm	5.0 mm	22	3 ml	1.5 ml	1 ml	0.25 ml
Code	Description	Quantity per box		Quantity per case		
218C	non sterile	500 pieces		3000 (6 boxes of 500 pcs.)		

NARROW STEM PIPETS



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw
83 mm	2.5 mm	22	4 ml	3.5 ml
Code	Description	Quantity per box		
203C	non sterile	1000 pieces		



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw	
150 mm	2.5 mm	28	4 ml	3.5 ml	
Code	Description	Quantity per box		Quantity per case	
202C	non sterile	500 pieces		5000 (10 boxes of 500 pcs.)	
202CS01	sterile individually wrapped	500 pieces			
202CS20	sterile 20 pieces per pack	25 packs of 20 (500 pcs total)			



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw	
63 mm	2.5 mm	25	1.2 ml	0.9 ml	
Code	Description	Quantity per box		Quantity per case	
209C	non sterile	500 pieces		3000 (6 boxes of 500 pcs.)	
209CS20	sterile 20 pieces per pack	25 packs of 20 (500 pcs total)			

BLOOD BANK PIPET



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw	Graduated up to	Grad. intvl.
155 mm	6.7 mm	22	5 ml	2.4 ml	2 ml	0.5 ml
Code	Description	Quantity per box		Quantity per case		
212C	non sterile	500 pieces		5000 (10 boxes of 500 pcs.)		
212CS01	sterile individually wrapped	500 pieces				
212CS20	sterile 20 pieces per pack	25 packs of 20 (500 pcs total)				

FINE TIP PIPETS



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw
104 mm	3.0 mm	50	1.5 ml	1.0 ml

Code	Description	Quantity per box	Quantity per case
222C	non sterile	400 pieces	2400 (6 boxes of 400 pcs.)
222CS01	sterile individually wrapped	500 pieces	
222CS10	sterile 10 pieces per pack	50 packs of 10 (500 pcs total)	
222CS20	sterile 20 pieces per pack	25 packs of 20 (500 pcs total)	



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw
150 mm	6.3 mm	66	7 ml	6 ml

Code	Description	Quantity per box	Quantity per case
223C	non sterile	400 pieces	4000 (10 boxes of 400 pcs.)
223CS01	sterile individually wrapped	500 pieces	
223CS20	sterile 20 pieces per pack	25 packs of 20 (500 pcs total)	



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw
144 mm	5.0 mm	65	5 ml	3.5 ml

Code	Description	Quantity per box	Quantity per case
224C	non sterile	500 pieces	5000 (10 boxes of 500 pcs.)
224CS01	sterile individually wrapped	500 pieces	
224CS20	sterile 20 pieces per pack	25 packs of 20 (500 pcs total)	

PADDLE PIPET



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw
130 mm	3.8 mm	33	3.0 ml	0.2 ml

Code	Description	Quantity per box	Quantity per case
220C	non sterile	500 pieces	3000 (6 boxes of 500 pcs.)

GENERAL PURPOSE PIPETS



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw
87 mm	5.0 mm	25	1.5 ml	1.0 ml

Code	Description	Quantity per box	Quantity per case
205C	non sterile	1000 pieces	6000 (6 boxes of 1000 pcs.)



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw
123 mm	5.0 mm	28	1.5 ml	1.0 ml

Code	Description	Quantity per box	Quantity per case
211C	non sterile	600 pieces	6000 (10 boxes of 600 pcs.)



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw
138 mm	5.0 mm	20	3.0 ml	1.5 ml


Code	Description	Quantity per box	Quantity per case
219C	non sterile	500 pieces	5000 (10 boxes of 500 pcs.)

EXTRA LARGE PIPETS



Total length	Stem diameter	Drops/ml	Capacity	Bulb draw
170 mm	4.0 mm	18	14 ml	10 ml

Code	Description	Quantity per box
207C	non sterile	1200 pieces

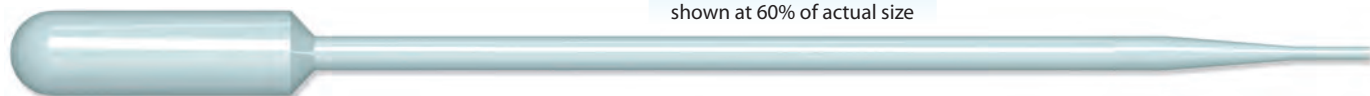

shown at 60% of actual size

Total length	Stem diameter	Drops/ml	Capacity	Bulb draw
300 mm	9.0 mm	22	23 ml	10 ml

Code	Description	Quantity per box	Quantity per case
227C	non sterile	100 pieces	600 (6 boxes of 100 pcs.)



NEW! DROPLESS PIPET


NO DROP COUNTING!

- Delivers accurate volume in one step
Simplifies rapid membrane test protocol
- Ideal for OTC kits
- Built-in procedural control for ***In Vitro*** diagnostic use e.g. home pregnancy test kits
- Choice of models to deliver different volumes: 40µl, 60µl, 80µl, 100µl and 120µl

**Just fill
the pipet
and release!**




pipet shown actual size

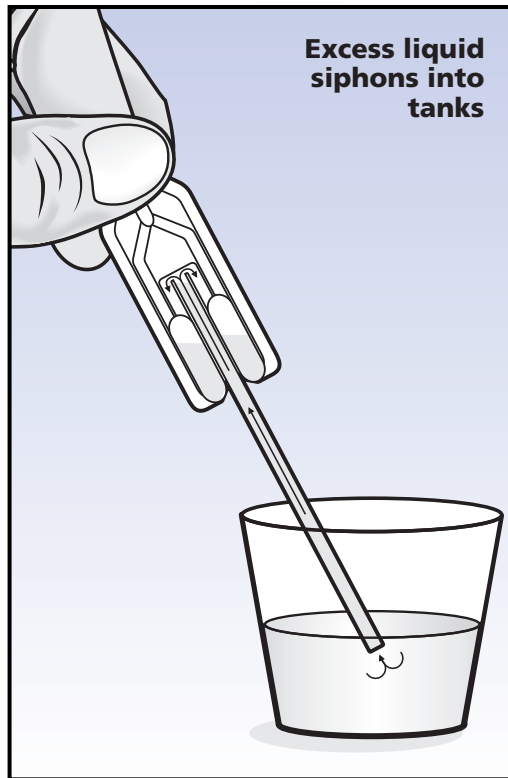

**Excess liquid
siphons into
tanks**

Squeeze and release bulb to draw up liquid. Excess liquid is siphoned automatically in two overflow tanks.

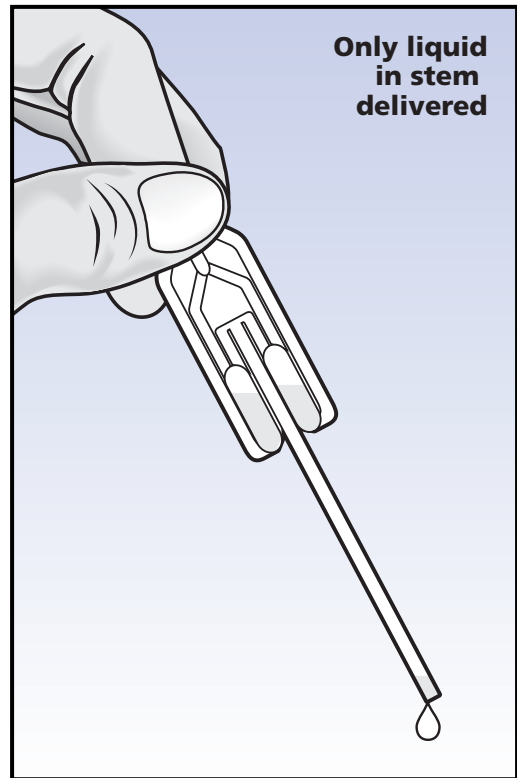

**Only liquid
in stem
delivered**

Squeeze bulb again to release ***precise volume*** just from the stem of the pipet.

Code	Description	Quantity per box	Quantity per case
214C	100µl pipet	2x500 pipets	6000 (6 boxes)



Copan Diagnostics Inc.
26055 Jefferson Ave. • Murrieta, California 92562 USA
Toll Free: (800)216-4016 (US & Canada)
Phone: (951)696-6957 Fax: (951)600-1832
E-Mail: info@copanusa.com

Local Distributor

www.copanusa.com